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declares that the object of the declaration is the product: Massage device

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conforms to essential requirement according to ECC Directive EMC 2014/30/UE, LVD 2014/35/UE, ROHS 2011/65/UE and DECRETO 15 October 2015, n. 206 applicable to this product. In addition, the object of the declaration described above, is in conformity with the relevant Union harmonization legislation:

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
This declaration of conformity is issued under the sole responsibility of the manufacturer.

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Additional information:

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Signed for and on behalf of:

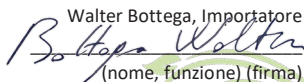Walter Bottega, Importatore

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Torino, 22/02/2018

ENGLISH

Dear customer,

thank you for choosing PressoEstetica® Mesis®. You will appreciate the high quality materials and features that set this model apart.

We hope that you find this machine comfortable to use and remind you to read the manual, the contraindications and warnings prior to use.

The product must be used in accordance with the instructions contained in this manual.

The manufacturer and distributor do not accept any liability for any use that differs from the one indicated in this manual.

PROFESSIONAL USE: professionals who purchase and use equipment for professional purposes are directly liable both legally and criminally. In particular, please note that professional use of equipment by aestheticians is regulated by Ministerial Decree-Law 1/90.

Laws and regulations, to which reference is made in this manual, are European and Italian.

VIMA s.r.l. Società a Socio Unico
Via Pigafetta 65, 10129 Torino
fax (+39) 011.0432817
email: info@mesis.eu
www.mesis.eu

Reproduction, even in part, of any photo or text contained within this manual is prohibited.

© Copyright VIMA s.r.l.

Index

Warnings24

Features PressoEstetica® Mesis®JoySense®31

PressoEstetica® Mesis® JoySense® 2.0 - mod: JOY2.032

PressoEstetica® Mesis® JoySense® 3.0 - mod: JOY3.033

PressoEstetica® Mesis® JoySense® instructions for Use34

Maintenance and Troubleshooting37

Warranty38

EU declaration of conformity for aesthetic use39

 **Warnings**

Before using the device, consult your doctor and read the instructions attentively for use and, in particular, safety warnings, thereby adhering to these.

Conserve the present manual for the entire service life for reference.

In the case of transferring the product to third parties, please also provide the documentation in its entirety.

After having removed the device from the packaging, check to ensure that the scope of supply is complete and check for any possible transport damage. If in doubt, do not use the device and contact the authorised service centre.

The packaging material is not a toy for children! Keep the plastic bags away from children; danger of suffocation!

Before connecting the device, check that the voltage data shown on the identification plate correspond to the specifications of electric mains available.

The identification label is located on the device.

The present device must only be used for the intended purpose it was designed for, or as a massage unit. Any other use is considered improper and therefore dangerous.

In order to ensure the best level of protection, it is advised to install a differential current device the rated operating differential current of which does not exceed 30mA (ideally 10mA). Ask your installer for advice concerning this.

Caution:

Do not use the present device near bathtubs, sinks or other receptacles containing water! Do not submerge the device into water!

The appliance can't be used by children aged 8 or less and by persons with reduced physical, sensory or mental abilities, or lacking experience or necessary knowledge, provided they are under surveillance or after they have received instructions regarding the safe use of the appliance and the understanding of the dangers inherent to it. Children must not play with the appliance. Cleaning and maintenance intended to be carried out by the user must not be carried out by children without supervision. Keep the appliance and its cord out of reach of children under 8 years.

DO NOT use the device with wet hands, damp or bare feet.

DO NOT pull the power cable or the device itself when disconnecting it from the power outlet.

DO NOT expose the device to moisture or the influences of weather (rain, sun).

Do not block the air openings of the appliance and do not place the appliance on a soft surface such as a bed or couch, where the openings may become blocked.

Keep the air openings free of any obstructions.

In the case of breakdown or malfunction of the device, turn it off and do not tamper with it.

Only contact an authorised technical service centre for any possible repairs.

If the power cord is damaged, it must be replaced by an authorised technical service centre in order to prevent any risk.

Use a moist cloth to clean the unit and check the contents of the present manual.

Please note that the apparatus must only be used indoors and in accordance with the operating instructions. If the power cord is faulty, it must be replaced by a VIMA s.r.l. customer service centre.

The electrical unit to which the device is connected has to be in conformity with the Ministerial Decree No. 37/08 (Italian Rules, then check the applicable regulations in your country) and, in particular, all the safety measures have to be in place for the unit - thermo-magnetic switches on the main power lines.

Please store and handle the product with care. It is prohibited to use the product if damage is present.

Do not position the device or its accessories near sources of heat.

Only use in domestic environments.

Photos and features the product refer to what is currently printed in the manual, differences of minor importance may be possible during the course of production. Possible differences shall have no influence on the safety or on the performance of the product.

ATTENTION: in order to avoid flame or shocks of various nature, do not expose the device to rain or liquid of any type.

The manufacturer shall provide information on the correct use of the items. Therefore, always adhere to these instructions that assure the correct and appropriate installation as well as functionality. For this reason, always keep these instructions for future reference.

Position the switch to 0 (zero) or OFF after discontinuing use.

All of the device's technical specifications must be borne in mind. Check the symbols that are on the ID tag of each item and the icons visible in the section (page 27-28-29-30) of the safety instructions.

This appliance may only be used indoors.

It is not suitable for installation in bathrooms.

The item is only suited to be directly connected to the mains using a proper power supply.

Electrical socket: in the case of damage, the electric cable must only be replaced by the manufacturer, the distributor or by an expert technician in order to prevent risks.

During installation ensure you do not make holes near the existing electrical wiring in the wall or ceiling.

Connecting cables must never be twisted and squashed into the space between the device and the supporting surface.

This product isn't suitable for children.

Do not treat:

- subjects with capillary fragility, arterial hypertension or heart disease;
- subjects with a pacemaker or an internal defibrillator.
- carriers of active implants such as pacemakers, insulin pumps or other devices;
- carriers of metallic prostheses;
- subjects with inflammatory processes in place, with skin lesions, with neoplasms;
- pregnant women;
- do not treat subjects with sensitive skin, who may be restricted to powers of very low disbursement;
- apply only on undamaged skin;
- never use the product directly on swollen or inflamed areas or skin rashes.

Consult with your medical specialist before use.

How to operate and further warnings:

The application time varies depending on how the massage is carried out and, as a rule, this can vary between 15 and 30 minutes.

Other than those listed above, attentively follow the indications, warnings and cautions for use specified within the present manual.

Mesis® JoySense machine is designed, manufactured and tested in compliance with the most stringent European safety standard regulations.

Any damage caused by using the device in exceptional circumstances (use outdoors, in wet areas, wet zones and exposure to direct sunlight ...) is not covered by the warranty terms.

The period of warranty shall begin from the date of purchase and shall only be deemed valid if a document of purchase, invoice or receipt is presented. The warranty shall lapse if the item has not been installed according to the specified instructions or if it is repaired or tampered with.

Furthermore, the manufacturer shall not assume any liability for damages caused due to misuse or incorrect application of the instructions.

We shall reserve the right to perform technical modifications or modifications to the design.

Environmental warnings

 In compliance to Legislative Decree No. 49 of March 14, 2014 (Italian regulatory reference).
„Implementation of Directive 2012/19/EU on waste electrical and electronic equipment (WEEE)“

The crossed bin symbol on the appliance indicates that the product at the end of its useful life must be collected separately from other waste. The user must therefore confer the equipment complete with essential components at the end of their life to the appropriate collection centers for electronic and electrotechnical waste, or return it to the seller when new equipment of an equivalent type is purchased, in one by one, or one zero for equipment with major side lower than 25cm. The adequate separate collection for the subsequent start of the equipment discharged to recycling, treatment and compatible environmental disposal helps to avoid possible negative effects on the environment and health and favors the recycling of the materials that make up the equipment. The illegal disposal of the product by the user involves the application of administrative sanctions in compliance to Legislative Decree No. 49 of March 14, 2014 (Italian regulatory reference). It is advisable to check the regulations in force in your country regarding the disposal of waste and Environmental warnings.

Cleaning

Always disconnect the electric power plug before cleaning or maintenance and in the case of non-use of the device.

Do not use stringent or abrasive detergents for cleaning the surface: use a damp cloth for cleaning the device when necessary.

If dirt is difficult to remove, use a light detergent and delicately rub the area in question.

The parts that come into contact with the skin must be clean between massages etc. using normal cleaning methods: gentle disinfectant on a soft cloth.

Disconnect from the power supply before cleaning.

It is recommended to properly sterilise and/or disinfect all the parts that come into contact with the

individual undergoing massage.

Clean items with a damp cloths and do not use solvents or abrasive products.

To clean the appliance, use a damp cloth and check the contents of this manual.

Avoid spreading liquids on all electrical parts.

Warnings

- Before use, carefully read the manual in its entirety;
- Keep this manual close to the device device, so you can consult it at any time;
- This manual explains in detail how to properly use the device in order to prevent danger and damage to persons or property.

Meaning of the symbols in this manual


 Warning	 Precautions
Failure to follow the instructions can result in serious injury or loss of life.	Failure to follow the instructions can result in the user causing injury to him/herself or damage to the device.
Indications	
 Indicates a restriction	 Is an indication to be followed properly
 Prohibits disassembly	 Indicates to unplug the power cord from the wall outlet

Warning

Do not use the product in the following cases:

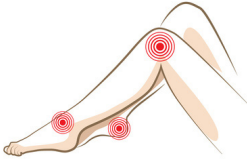People suffering from both internal and external inflammation or wounds at risk of bleeding must not use the product including those suffering from purulent inflammation.

Do not use in case of general malaise, flu or illness.

Contraindications

- Pulmonary oedema
- Thrombophlebitis
- Heart failure
- Infections
- Deep vein thrombosis
- Pulmonary embolism episodes
- Wounds, lesions or tumours in the area of application or in proximity to areas where increased venous and lymphatic return could be harmful
- Carriers of active implants such as pacemakers, insulin pumps or other devices;
- Carriers of metallic prostheses;
- Subjects with inflammatory processes in place, with skin lesions, with neoplasms;
- Pregnant women.

Do not treat:

- carriers of active implants such as: pacemakers, insulin pumps or other devices;
- carriers of metallic prostheses;
- subjects with inflammatory processes in place, with skin lesions, with neoplasms;
- pregnant women;
- do not treat subjects with sensitive skin, who may be restricted to powers of very low disbursement;
- apply only on undamaged skin.

Main contraindications

Pulmonary oedema. Thrombophlebitis. Heart failure. Infections. Deep vein thrombosis. Pulmonary embolism episodes. Wounds, lesions or tumours in the area of application or in proximity to areas where increased venous and lymphatic return could be harmful. Carriers of pacemakers. Lymphangitis. Plexopathy and neuropathy.

Consult your doctor to check whether you have any disease which is incompatible with the Mesis® JoySense massage. A non-exhaustive list of incompatible diseases are listed in the sections „Contraindications“, „Do not treat“, and „Main contraindications“ just mentioned.

Warnings

Warning

	<p>If you feel unwell or experience pain during use, do not use the device. Discontinue use immediately.</p>		<p>Do not use the device for more than 15 minutes during the first massage. Gradually increase the time in subsequent sessions depending on how well you tolerate the massage. Begin the massage starting with minimal pressure.</p>
	<p>Use for children is prohibited. Keep out of the reach of children. Risk of electric shock or other damage.</p>		<p>Do not use water or other substances that could come into contact with internal parts of the device.</p>
	<p>In case of failure or malfunction, do not open or attempt to repair the device. Risk of electric shocks and injuries to persons or property.</p>		<p>Your hands must be dry when connecting or disconnecting the power cord from the power supply (socket). Risk of electric shocks and injuries to persons or property. Risk of fire.</p>
	<p>Use only single-phase 220V current. After using the device, ensure you disconnect the power cord from the power supply properly. Risk of electric shocks and injuries to persons or property. Risk of fire.</p>		<p>Leave the power cord free, do not place it on top of, squash it against or attach it to any object. Risk of shock or electrical leakage. Risk of fire.</p>
	<p>Use a plug socket for the exclusive use of the device. Too many plugs could cause the start of a fire, electrical leakage and electric shocks.</p>		<p>Do not hit the device or let it fall. This could cause serious damage to the device.</p>
	<p>Do not insert any object or product into the device. Risk of damage to the device and to people.</p>		<p>Do not place the device and its accessories next to sources of heat. It could cause serious damage to the device and/or accessories</p>

Warnings

Warning

	<p>When wearing the leggings, abdominal or armband, do not move or walk.</p>		<p>Do not place the device and its accessories in direct sunlight or near the stove (gas, electric or any other fuel). This may cause deterioration, deformation and discolouration of the product.</p>
	<p>Do not use the device for a purpose other than massage. Risk of damage to people and the device.</p>		<p>Do not bend, knot, cut or break the power cord.</p>
	<p>Do not place any object on top of the device. Do not place weights on the machine or on any of its connections, nor on electric cables and connecting air tubes. It may cause changes in the pressure.</p>		<p>If the machine has not been used for a while or is used in a cold place, it may not work properly, therefore, once switched on, leave to work at a low pressure for at least 3 inflation cycles. Do not leave the machine in a damp place or in the rain.</p>
	<p>Do not use water, oil, gasoline or chemicals to clean the accessories (leggings, armbands, etc.). Use a damp cloth. It could damage them, shortening their lifespan.</p>		<p>Do not use the accessories directly on the skin. Wear suitable undergarments underneath the applicators.</p>
	<p>Do not use water, oil, gasoline or chemicals to clean the accessories. Use a damp cloth.</p>		<p>Do not bend, knot, or cut the air tubes.</p>

Features PressoEstetica® Mesis®JoySense®

- The PressoEstetica® Mesis® JoySense® are **sequential** and intermittent;
- **Program** - chamber inflation mode;
- **Chambers or Sectors** - elements that comprise the applicator and that fill with air, in the PressoEstetica® Mesis®JoySense® there are 5 chambers for every applicator;
- **Cycle** – one cycle is equivalent to the inflation of all chambers, from the first to the last regardless of the program used. After the last chamber has inflated, there is a pause and then a new cycle starts;
- **Applicator** – pneumo relaxation element composed of several chambers that inflate sequentially regulated by the chosen program (leggings, armband and abdominal/buttocks band);
- **Time** - duration of the session (adjustable by 15 – 25 - 30 minutes);
- **Pressure** - the quantity of air blown in the chamber, the pressure is adjustable and continues on each chamber.
- **Connector** - connecting tubes between the machine and the applicators. Thanks to "one touch connect", inserting the connecting tubes for the air to the machine is easy and quick.

The massage is **sequential**, the pressure starts to inflate the chamber located at the end of the limb (e.g. in the legs it starts from the feet) and once the chambers are inflated, reaches the base of said limb (e.g. In the legs to the groin).

PressoEstetica® Mesis® JoySense® 2.0 - mod: JOY2.0

Main Technical Features

- **2 programs, 5 CHAMBERS**
- Adjustable massage time
- Input: 220 - 240 Vac, 50/60Hz
- Power: 30 W
- Size 340X290X300mm
- Device weight: 1,8Kg
- 24 month warranty
- EU declaration of conformity for aesthetic use

Importatore: VIMA s.r.l.
 Via Pigafetta 65, 10129 Torino, Italia
 info@mesis.eu - www.mesis.eu - Fax (+39) 011.0432817
 Modello: JoySense 2.0 (JOY2.0)
 Nome commerciale: JoySense
 S/N
 Input: 220 - 240 Vac, 50/60Hz
 Power: 30 W
 Fabbricato in Korea


START-TIMER:
 Press to start the massage
 Press after switching on to vary the time

PROGRAM SELECTION:
 Slim: all full
 Silhouette: sequential massage

SWITCHING ON THE MACHINE:
 Press to turn the machine on

ADJUSTING THE INTENSITY:
 Turn the control knob to increase or decrease the intensity

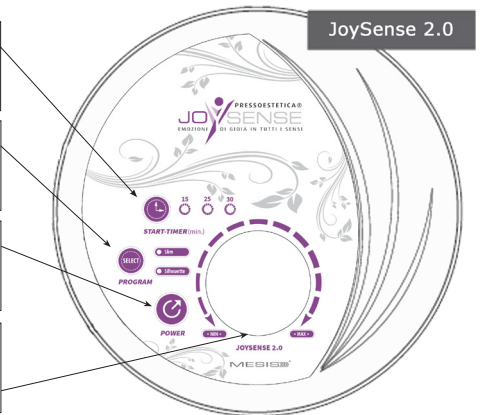Accessories, the equipment varies depending on the chosen configuration

PressoEstetica® Mesis® JoySense® 3.0 - mod: JOY3.0


Main Technical Features

- **3 programs, 5 CHAMBERS**
- Adjustable massage time
- Input: 220 - 240 Vac, 50/60Hz
- Power: 30 W
- Size 340X290X300mm
- Device weight: 1,8Kg
- 24 month warranty
- EU declaration of conformity for aesthetic use

Importatore: VIMA s.r.l.
 Via Pigafetta 65, 10129 Torino, Italia
 info@mesis.eu - www.mesis.eu - Fax (+39) 011.0432817
 Modello: JoySense 3.0 (JOY3.0)
 Nome commerciale: JoySense 3.0
 S/N
 Input: 220 - 240 Vac, 50/60Hz
 Power: 30 W
 Fabbricato in Korea

- START-TIMER:**
 Press to start the massage
 Press after switching on to vary the time
- PROGRAM SELECTION:**
 Slim: all full
 Silhouette: sequential massage
 Beauty: deep sequential massage
- SWITCHING ON THE MACHINE:**
 Press to turn the machine on
- ADJUSTING THE INTENSITY:**
 Turn the control knob to increase or decrease the intensity

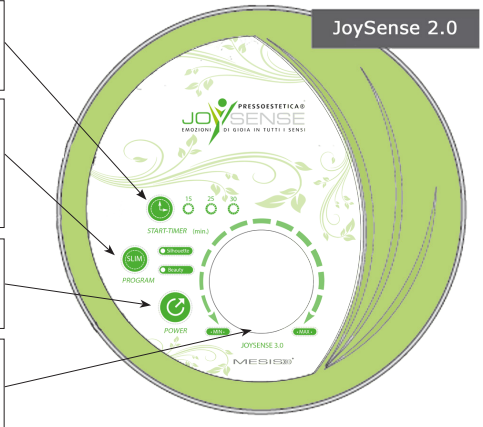Accessories, the equipment varies depending on the chosen configuration

Machine

Air hose plug-in area

Acupressure soles

Leg cuffs

Addominal cuff

Air connettor hoses


Extenders

Arm cuff


PressoEstetica® Mesis® JoySense® instructions for Use

1. Connect the power cord.

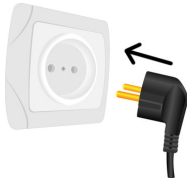
2. Connect the applicator's connector to use to the machine

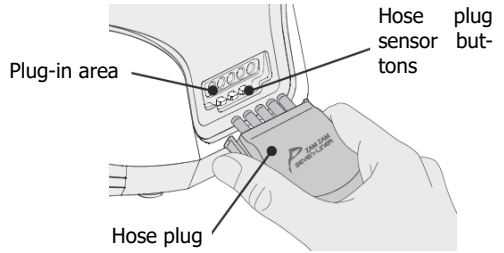3. Connect the applicators to the related connectors.

Leg cuffs



- connect the attachments of the the double connector to the leg cuffs.

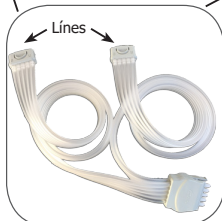Double connetor

Leg cuffs and abdominal/buttocks belt

- if you use the abdominal belt and leg cuffs, connect the double connector to the machine after that, you will have to connect one line of double connector to the splitter and the other one, to abdominal belt. Finally, connect the two lines of the splitter to the leggings.

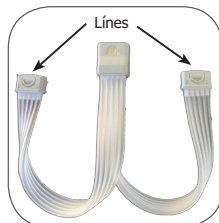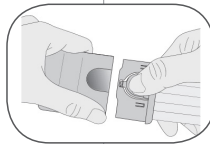Splitter

Abdominal/buttocks belt

- if you use only the abdominal belt, connect the attachments of the single connector to the belt.

*Arm cuff*

- connect the attachments of the single connector to the cuff.

4. Put on the applicators that you wish to use.

Leg cuffs

- Wear the leggings and close them with the zip, place the leggings with the attachments turned inwards.
- Inflation must start at the feet and gradually reach the groin.

Abdominal/buttocks belt

- Wear the abdominal belt on the buttocks-hips-abdomen and fasten it with the Velcro.
- Inflation must start at the top part and end at the base of the abdominal belt.

Arm cuff

- Wear the armband on and fasten it with the zip.
- Inflation must start from the hand and move sequentially to the shoulder.

5. Press „POWER“ to turn the machine on.

The lamp is turned on the last used program.


6. Press „PROGRAM“ to select the program.


In the model JoySense® 3.0, when you select the program „SLIM“ remain lit the lights of both programs „SILHOUETTE“ and „BEAUTY“.

Programs JoySense® 2.0

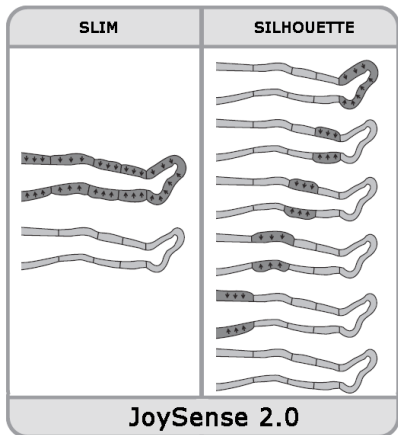Programs JoySense® 3.0

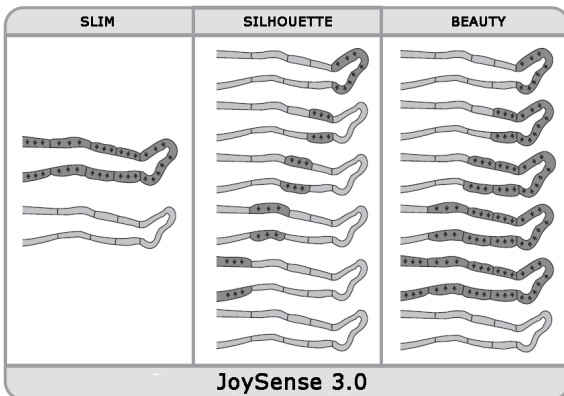
7. Press „START/TIME“ to set the duration of massage. Press it again to change the preset time (15 - 25 - 30 mins) .

8. Turn the knob to adjust the intensity.

- To increase the intensity, turn to the right
- To decrease the intensity rotate to the left

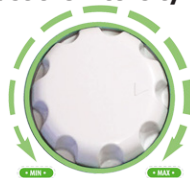N.B. When you decrease the intensity during the „SLIM“ program, inflation is stopped for about 5 seconds. In the other programs „SILHOUETTE“ and „BEAUTY“, when you decrease the intensity, the compressor automatically sets a shorter time inflation of the rooms.

Inflation starts only after have rotated the knob of the intensity adjustment.

Maintenance and Troubleshooting

The accessories cannot be repaired, as such take care of them.

Do not place the leggings, armband, abdominal belt or other accessories close to direct or indirect sources of heat.	
Do not place the equipment in damp or wet areas. Do not soak the equipment in any way, do not pour any liquid over the equipment.	
Do not wash the equipment and accessories with water or detergents. Simply wipe with a damp cloth.	
Do not place the equipment and its accessories in very cold areas or areas below zero degrees centigrade. It could cause serious damage.	
If you don't use the device for a length of time, place it back in its box.	

Environment of use and safekeeping

- **Use and maintain the device and its accessories in a clean environment and at a temperature between 10 – 40 °C, and humidity of 30 – 85%;**
- **Do not use outdoors.**

Problem found	Possible solution
The machine does not turn on	<ul style="list-style-type: none"> • Check that it is correctly connected to the mains; • Check that you have inserted the air tube connector to the machine; • Check the condition of the power cord. It must be intact.
If you hear a strange hissing sound similar to an air outlet	<ul style="list-style-type: none"> • Check that the tubes, accessories and connectors are intact; • Check that all the connections (to the machine and the applicators) have been correctly inserted; • Check that the air connection tubes are not bent, squashed or crushed.
The air does not flow to the applicators, there is insufficient intensity	<ul style="list-style-type: none"> • Check that the connectors to the machine and the applicators are correctly connected; • Check that the air connection tubes are not bent, squashed or crushed; • Check that the chambers are inflated progressively and according to the selected program; • Make sure you set the pressure control correctly; • Verify that the actual air outlet from the unit is present.
The air continues to flow into the chambers and you hear an air leak in the applicators	<ul style="list-style-type: none"> • Check for holes or damage in the applicators' air chambers. If no anomalies are detected, switch off the device, switch it back on and check that it is working properly again

 **Warranty**

- The device is guaranteed for the first user for a period of twenty four (24) months, three (3) months for the applicators and accessories subject to wear from the date of purchase against defects in materials or manufacture, provided that it is used properly and maintained in good working condition.
- The warranty is valid, if the user uses the device for professional purposes, for twelve (12) months for the machine and three (3) months for the applicators and accessories subject to wear, from the date of purchase against defects in materials or manufacture, provided that it is used properly and maintained in good working condition.

The warranty's validity is limited in the following cases:

- No warranty for accessories and materials to be deemed as "consumables"
- This Warranty is valid and effective in the country in which the Product was purchased. In the event that the product was purchased in an EU country, the warranty is still valid in all EU countries.
- Problems due to dropping the device or to improper use of the voltage.
- If the device has been opened or modified by unauthorised personnel.
- If the device has been used in unsuitable environments.
- If the device has not been used in accordance with the instructions contained in this manual.
- Warranty on the device: 2 years.
- Warranty on accessories and applicators: 3 months.

To obtain the warranty service, the user must observe the following warranty clauses:

- The products must be delivered for repair in the original packaging.
- The product warranty will be subject to the production of a tax document (sales receipt, payment receipt, sales invoice or mail-order sales document via an online shop or catalogue), attesting the product's date of purchase.
- Repair will not be performed on the original date of the warranty's expiry and it will not constitute any renewal or extension of the same.
- If, at the time of repair, there is no fault in the product found, expenses related to the technical verification shall be charged.
- The warranty shall lapse if malfunction has been determined to be due to impact, falls, incorrect or improper use of the product, use of a power supply or charger that is not the original, accidents, alterations to the product, substitution/removal of the warranty seals and/or tampering with the product.
- Moreover, the same warranty shall not cover damage caused during transport through use of unsuitable packaging.
- The warranty does not cover the impossibility to use the product, other incidental or consequential costs or other expenses incurred by the purchaser.

PRODUCT SUPPORT

In the event of failure or malfunction of the device, contact your dealer.

Distributed by:

VIMA s.r.l. Società a Socio Unico
Via Pigafetta 65, 10129 Torino, Italy
P.I. 08218930017
email: info@mesis.eu

EU declaration of conformity for aesthetic use

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Da: Welbutech Co., Ltd.