**MOPEDIA®**

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---

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**INDICE**

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**MOPEDIA®**

MA XMP370-MP380-MP381-MP390-MP391 02\_A\_ENG\_09 2020

**HOSPIYAL BEDS THREE-JOINTED ELECTRIC BEDS**

---

**INSTRUCTION MANUAL**

ENGLISH

## INDEX

<b>1. CODES.....</b>	<b>PAG.3</b>
<b>2. INTRODUCTION.....</b>	<b>PAG.3</b>
<b>3. INTENDED USE.....</b>	<b>PAG.3</b>
<b>4. DECLARATION OF CONFORMITY .....</b>	<b>PAG.3</b>
4.1 Reference standards and directives.....	pag. 4
<b>5. GENERAL WARNINGS .....</b>	<b>PAG.4</b>
<b>6. SYMBOLS .....</b>	<b>PAG.4</b>
<b>7. GENERAL DESCRIPTION.....</b>	<b>PAG.5</b>
7.1 Bed structure MP370 (three-jointed electric hospital bed mod. GERANIO).....	pag. 5
7.2 Electrical kit for MP370 (three-jointed electric hospital bed mod. GERANIO) .....	pag. 5
7.3 Bed structure MP380 (three-jointed electric lifting hospital bed mod. GERBERA).....	pag. 5
7.4 Electrical kit for MP380 (three-jointed electric lifting hospital bed mod. GERBERA).....	pag. 6
7.5 Bed structure MP381.....	pag. 6
7.6 Electrical kit MP381.....	pag. 6
7.7 Bed structure MP390 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA).....	pag. 6
7.8 Electrical kit for MP390 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA).....	pag. 7
7.9 Bed structure MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA).....	pag. 7
7.10 Electrical kit for MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA).....	pag. 7
<b>8. ASSEMBLY .....</b>	<b>PAG.8</b>
8.1 MP370 (three-jointed electric hospital bed mod. GERANIO) .....	pag. 8
8.2 MP380 and MP381 (three-jointed electric lifting hospital bed mod. GERBERA).....	pag. 11
8.3 MP390 and MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA) .....	pag. 14
8.4 Assembly of the handles on the sections (for MP370/MP380/MP381/MP390/MP391) .....	pag. 16
<b>9. BEFORE EVERY USE .....</b>	<b>PAG.16</b>
<b>10. WARNINGS ON USE.....</b>	<b>PAG.17</b>
<b>11. HOW TO USE.....</b>	<b>PAG.17</b>
11.1 MP370 (three-jointed electric hospital bed mod. GERANIO) .....	pag. 17
11.2 MP380-MP381 (three-jointed electric lifting hospital bed mod. GERBERA) .....	pag. 17
11.3 MP390-MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA) .....	pag. 18
<b>12. MAINTENANCE .....</b>	<b>PAG.19</b>
<b>13. CLEANING AND DISINFECTION .....</b>	<b>PAG.19</b>
13.1 Cleaning.....	pag. 19
13.2 Disinfection.....	pag. 19
<b>14. CONDITIONS OF DISPOSAL .....</b>	<b>PAG.19</b>
14.1 General conditions of disposal .....	pag. 19
14.2 INSTRUCTIONS FOR CORRECT DISPOSAL IN ACCORDANCE WITH EUROPEAN DIRECTIVE 2002/96/EC:.....	pag. 19
<b>15. DECLARATION OF ELECTROMAGNETIC COMPATIBILITY .....</b>	<b>PAG.19</b>
<b>16. REPLACEMENT PARTS AND ACCESSORIES .....</b>	<b>PAG.22</b>
<b>17. TROUBLESHOOTING.....</b>	<b>PAG.23</b>
<b>18. TECHNICAL SPECIFICATIONS.....</b>	<b>PAG.24</b>
18.1 Technical specifications for all models .....	pag. 24
18.2 Dimensions .....	pag. 24
<b>19. WARRANTY .....</b>	<b>PAG.26</b>
<b>20. REPAIRS .....</b>	<b>PAG.26</b>
<b>21. NON-LIABILITY CLAUSE .....</b>	<b>PAG.26</b>

## CE | Class Medical Device

### REGULATION (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 5 April 2017 concerning medical devices

#### 1. CODES

<b>MP370</b>	three-jointed electric hospital bed, fixed height
<b>MP380</b>	three-jointed electric hospital bed, variable height
<b>MP381</b>	three-jointed electric hospital bed with slats, variable height
<b>MP390</b>	three-jointed electric "Trendelenburg" hospital bed, variable height
<b>MP391</b>	three-jointed electric, easy "Trendelenburg" hospital bed, variable height

#### 2. INTRODUCTION

Thank you for choosing the **three-jointed electric hospital bed of MORETTI S.p.A.** Their design and quality guarantees unbeatable comfort, safety and reliability. **The three-jointed electric hospital beds in the MOPEDIA range by MORETTI S.p.A.** have been designed and manufactured to satisfy all your needs for practical and correct use. This user manual contains some suggestions on how to use the wheelchair properly and safely. Please read this manual in full before you use the device. Should you have any queries, please contact your dealer for suitable advice and assistance.

#### 3. INTENDED USE

MOPEDIA hospital beds are intended for patients admitted for health reasons at hospitals and/or healthcare structures and entities or for use at home.

#### WARNING

- Do not use the product for any purpose other than that specified in this manual.
- Moretti S.p.A. declines any and all liability for damages resulting from improper use of the device or any use other than that specified in this manual.
- The manufacturer reserves the right to modify the device and its relevant manual without prior notice, in order to improve the product's features.

#### 4. DECLARATION OF CONFORMITY

Moretti S.p.A. declares under its sole responsibility that the product made and traded by Moretti S.p.A. and belonging to the group of MEDICAL BEDS - MOPEDIA complies with the provisions of the regulation 2017/745 on MEDICAL DEVICES of 5 April 2017.

For this purpose, Moretti S.p.A. guarantees and declares under its sole responsibility what follows:

1. The devices satisfy the requirements of general safety and performance requested by the Annex I of regulation 2017/745 as laid down by the Annex IV of the above mentioned regulation.
2. The devices ARE NOT MEASURING INSTRUMENTS.
3. The devices ARE NOT MADE FOR CLINICAL TESTS.
4. The devices are packed in NON-STERILE BOX.
5. The devices belong to class I in accordance with the provisions of Annex VIII of the above mentioned regulation
6. Moretti S.p.A. provides to the Competent Authorities the technical documentation to prove the conformity to the 2017/745 regulation, for at least 10 years from the last lot production.

**Note:** Complete product codes, the manufacturer registration code (SRN), the UDI-DI code and any references to used regulations are included in the EU declaration of conformity that Moretti S.p.A. releases and makes available through its channels.

#### 4.1 Reference standards and directives

To ensure safety standards for users of the product, Moretti S.P.A tests and repeats the legislation: UNI CEI 60601-2-52.

### 5. GENERAL WARNINGS

#### WARNING

Take extreme care in the presence of moving parts that could trap limbs and cause injury.

- Refer to this manual for correct use of the device, paying particular attention to maintenance and to the precautions before every use.
- Consult your doctor or therapist for correct use of the device.
- Keep the packaged product away from sources of heat, as the packaging is cardboard.
- The lifetime of the device depends on wear and tear of non-repairable and/or non-replaceable parts.
- Always take care when children are present.
- Do not exceed the maximum load for safe use.
- Never clean with aggressive products such as: Ammonia, bleach, acids and abrasive detergents.
- The user and/or the patient will have to report any serious accident that have occurred ..... related the device to the manufacturer and appropriate authority of the State which the .... user and/or patient belongs to.

### 6. SYMBOLS

	Product code
	Unique Device Identification
	CE mark
	Manufacturer
	Batch Lot
	Read the instruction manual
	Medical Device
	Conditions of disposal
	Attention, should read the instructions
	Production dates
	Waste Electrical & Electronic Equipment (WEEE)

**IP54** Protection against dust and water splashes

B Type

Insulation class

Maximum permissible load

Maximum Load

## 7. GENERAL DESCRIPTION

The hospital beds are formed of the following parts:

### 7.1 Bed structure MP370 (three-jointed electric hospital bed mod. GERANIO)

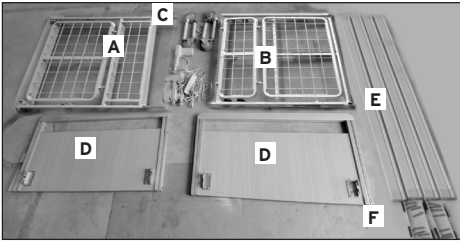
Fig. 1

#### List of parts

- A N°1 mesh frame head raiser part
- B N°1 mesh frame leg raiser part
- C N°4 supports with castors with brake
- D N°2 headboard/footboard
- E N°4 side slats
- F N°4 side slides

### 7.2 Electrical kit for MP370 (three-jointed electric hospital bed mod. GERANIO)

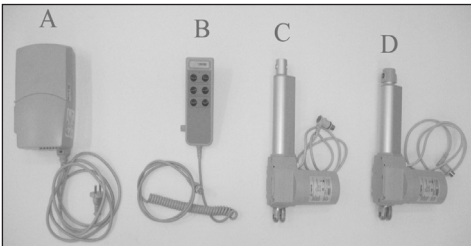
Fig. 2

#### List of parts

- A N°1 electronic control unit
- B N°1 remote control
- C N°1 head raiser electric actuator
- D N°1 leg raiser electric actuator

### 7.3 Bed structure MP380 (three-jointed electric lifting hospital bed mod. GERBERA)


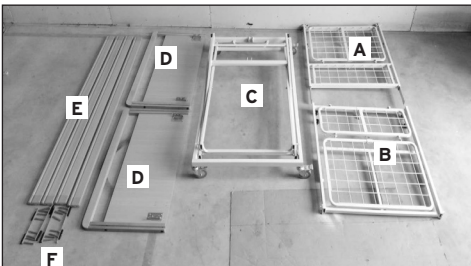
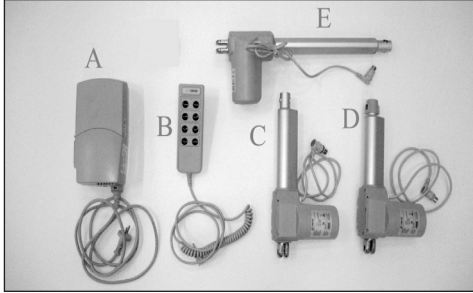
Fig. 3

#### List of parts


- A N°1 mesh frame head raiser part
- B N°1 mesh frame leg raiser part
- C N°1 lifting frame with wheels
- D N°2 headboard/footboard
- E N°4 side slats
- F N°4 side slides

**7.4 Electrical kit for MP380 (three-jointed electric lifting hospital bed mod. GERBERA)**

**List of parts**

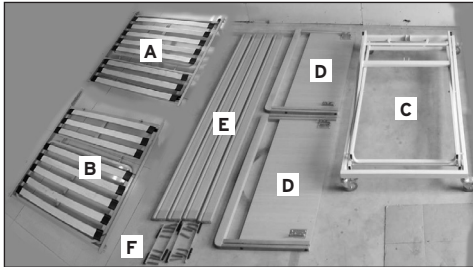
- A N°1 electronic control unit
- B N°1 remote control
- C N°1 head raiser electric actuator
- D N°1 leg raiser electric actuator
- E N°1 lifting electric actuator

**Fig. 4**
**7.5 Bed structure MP381**

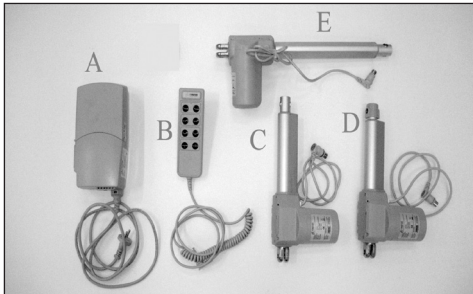**List of parts**

- A N°1 mesh frame head raiser part
- B N°1 mesh frame leg raiser part
- C N°1 lifting frame with wheels
- D N°2 headboard/footboard
- E N°4 side slats
- F N°4 side slides

**Fig. 5**
**7.6 Electrical kit MP381**

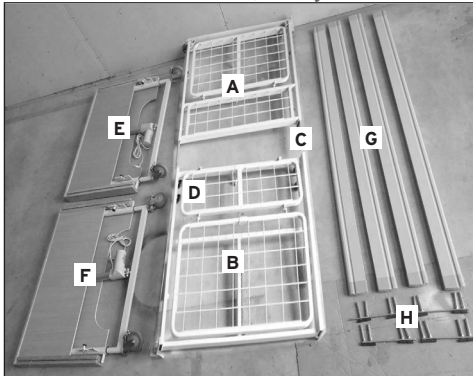**List of parts**


- A N°1 electronic control unit
- B N°1 remote control
- C N°1 head raiser electric actuator
- D N°1 leg raiser electric actuator
- E N°1 lifting electric actuator

**Fig. 6**
**7.7 Bed structure MP390 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA)**

**List of parts**

- A N°1 mesh frame head raiser part
- B N°1 mesh frame leg raiser part
- C N°2 frame connection rectangular tubes
- D N°4 knobs to lock the frame
- E N°1 headboard
- F N°1 footboard
- G N°4 wooden slats for the sides
- H N°4 plastic slides for the sides

**Fig. 7**

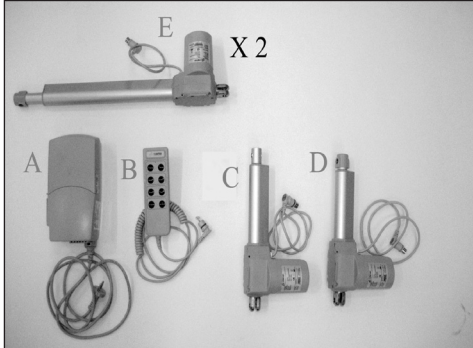**7.8 Electrical kit for MP390 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA)**


**Fig. 8**

**List of parts**

- A N°1 electronic control unit
- B N°1 remote control
- C N°1 head raiser electric actuator
- D N°1 leg raiser electric actuator
- E N°2 headboard/footboard electric actuators

**7.9 Bed structure MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA)**

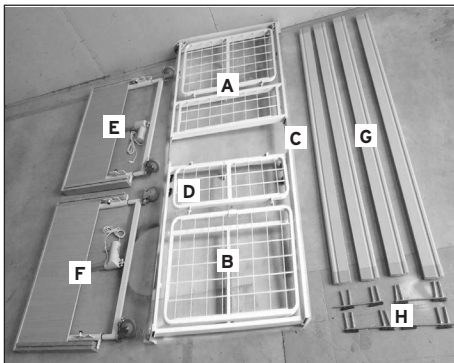
**Fig. 9**

**List of parts**

- A N°1 mesh frame head raiser part
- B N°1 mesh frame leg raiser part
- C N°2 frame connection rectangular tubes
- D N°4 knobs to lock the frame
- E N°1 headboard
- F N°1 footboard
- G N°4 wooden slats for the sides
- H N°4 plastic slides for the sides

**7.10 Electrical kit for MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA)**

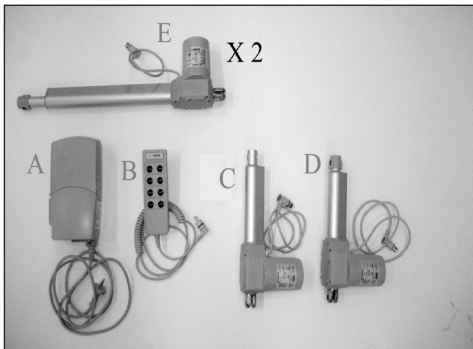**Fig. 10**

**List of parts**

- A N°1 electronic control unit
- B N°1 remote control
- C N°1 head raiser electric actuator
- D N°1 leg raiser electric actuator
- E N°2 headboard/footboard electric actuators

## 8. ASSEMBLY

### 8.1 MP370 (three-jointed electric hospital bed mod. GERANIO)

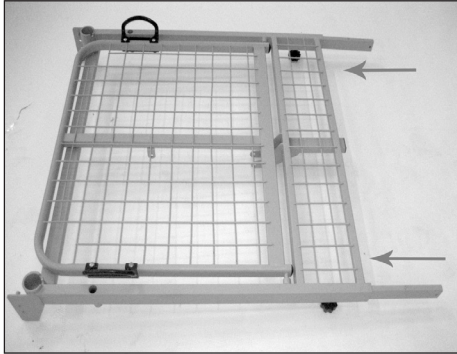
Fig. 11

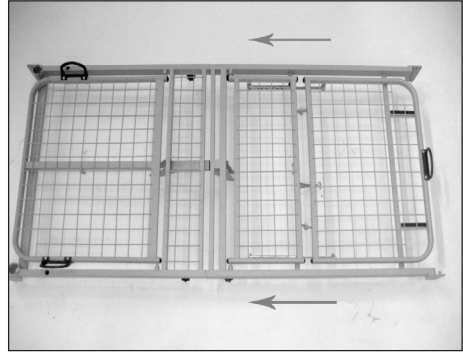
Fig. 12

1. Insert the two rectangular tubes into the frame at the head raiser end and push them right in (Fig. 11)
2. Join the two parts of the mesh at this point, aligning them and locking them in place with the four knobs (Fig. 12)
3. Assemble the four wheel supports at the respective corners and lock them in place with the bolt provided (Fig. 13)

Fig. 13

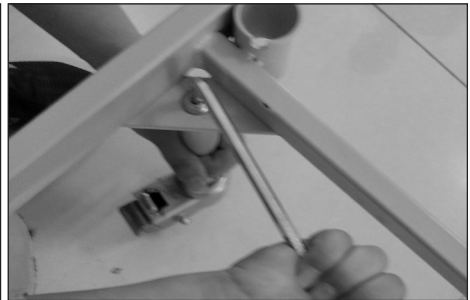
Fig. 14

4. Screw the wheels onto the supports (it is recommended to lock the brake on each wheel to perform this operation, to make it easier to screw them on by hand). Tighten completely (Fig. 14)

5. Insert the wooden headboard and footboard (two identical elements) in the supports on the frame, making sure that the automatic locking pin is correctly inserted into the hole (Fig. 15)

**Fig. 15**


6. Prepare the sides for assembly in the following manner:

- Take a pair of slats and position them parallel to the floor

N.B. THE SIDES HAVE A SMOOTH SIDE, A LOWER PART AND A ROUNDED SIDE THAT IS THE UPPER PART. POSITION THE SMOOTH SIDE OF BOTH IN THE SAME DIRECTION

**Fig. 16**

- Insert the plastic slides as a connection element of the sides on both sides, taking care to ensure that the pointed triangular part remains in the same direction as the rounded part of the side (Fig. 16)
- Carry out the same operation for the other two slats, in a manner to form the pair of sides

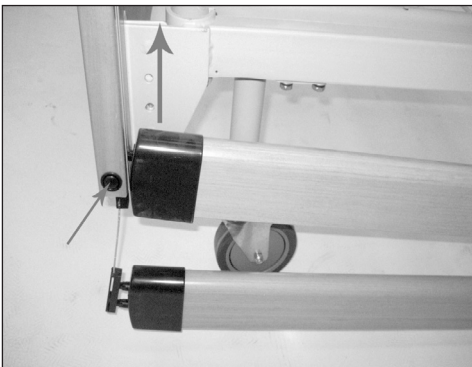7. Insert the two pairs of slats into the lower part of the guides on the wooden headboard and footboard. Perform this operation first on one side and then on the other, pressing the lock button to insert the slat in the guide, then releasing the button to lock them in place (Fig. 13)

**Fig. 17**

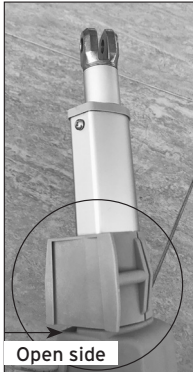
Fig. 18

8. Install the control unit on the actuator with the shortest cable. Perform the following operations before installing the control unit on the actuator:

- Connect the power cable to the control unit;
- Insert the plastic support onto the actuator, as shown in Fig. 18, with the open side of the tapered part facing downwards;
- Install the control unit on the actuator by sliding it into the guides and locking it in place using the plastic support installed in the previous point. Fig. 19;
- Install the actuator with the control unit on the two welded pins underneath the mesh structure at the headboard end. As follows:
  - Align the hole on the actuator with the one on the brackets of the bed section (Fig. 20);
  - Insert the quick-release pin into the hole (ensure the bracket is properly positioned inside the actuator) (Fig. 21);
  - Close the stop mechanism on the quick-release pin (Fig. 22).

Fig. 19

9. Assemble the other actuator underneath the structure of the mesh at the leg raiser end, repeating the operations described in the previous point.

Fig. 20

Fig. 21

Fig. 22

10. Make the electrical connections on the control unit, as follows:

- Remove the plastic cover from the control unit and insert a flat screwdriver into the small window on the cover (Fig. 23);
- connect the cable of the actuator installed on the control unit on the "back" channel;
- connect the cable of the actuator on the foot raiser side on the "foot" channel;

Fig. 23



- A - Not connected
- B - Not connected
- C - Leg raiser electric actuator
- D - Head raiser electric actuator
- E - Remote control

- Position the cables in the cable housings present on the control unit and reposition the plastic cover

At this point, the bed is assembled and ready for use.

**8.2 MP380 and MP381 (three-jointed electric lifting hospital bed mod. GERBERA)**

1. Insert the two rectangular tubes into the frame at the head raiser end and push them right in (Fig. 24)

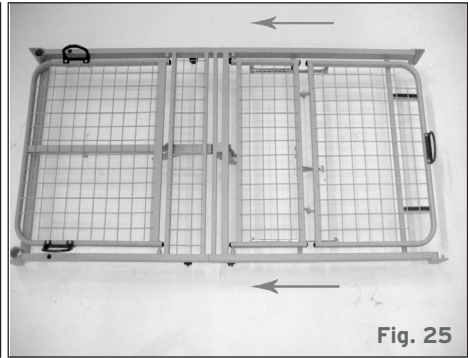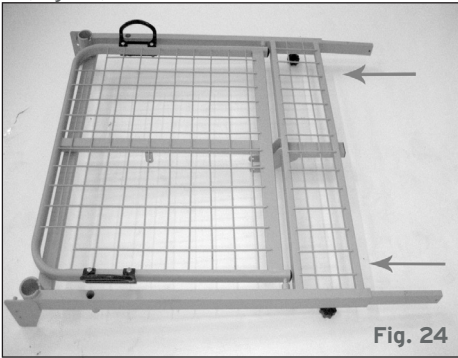Illustrative images

2. Join the two parts of the mesh at this point, aligning them and locking them in place with the four knobs (Fig. 25)

3. Connect the lifting frame to the mesh in the following manner:

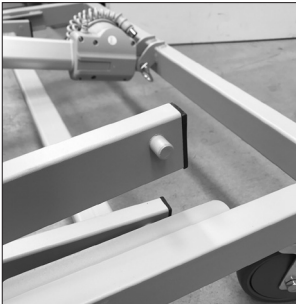
Fig. 26

Fig. 27

- Place the mesh frame and lifting frame in the correct mutual position. Insert the two black bearings into the two pins of the lifting frame, as shown in Fig. 26 and Fig. 27.

Fig. 28

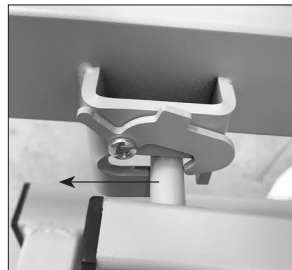
Fig. 29

- Assemble the frames by sliding the bearings inside the guides of the mesh frame (Fig. 28);
- Slide the mesh frame until the pins of the lifting frame are engaged with the hook on the mesh frame (Fig. 29).

4. Use the screws and wrenches provided to assemble the metal supports on the headboard and footboard. Assemble the supports in a manner that the open part of the hole is facing downwards.

5. Insert the wooden headboard and footboard (two identical elements) in the supports on the frame, making sure that the clasp rotates and locks the pin in place (Fig. 30).

Fig. 30

6. Prepare the sides for assembly in the following manner:

- Take a pair of slats and position them parallel to the floor

Fig. 31

N.B. THE SIDES HAVE A SMOOTH SIDE, A LOWER PART AND A ROUNDED SIDE THAT IS THE UPPER PART. POSITION THE SMOOTH SIDE OF BOTH IN THE SAME DIRECTION

- Insert the plastic slides as a connection element of the sides on both sides, taking care to ensure that the pointed triangular part remains in the same direction as the rounded part of the side (Fig. 31)
- Carry out the same operation for the other two slats, in a manner to form the pair of sides

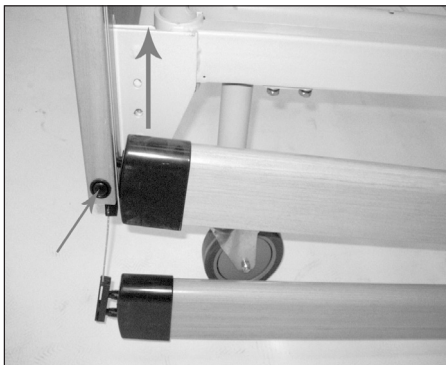
Fig. 32

7. Insert the two pairs of slats into the lower part of the guides on the wooden headboard and footboard. Perform this operation first on one side and then on the other, pressing the lock button to insert the slat in the guide, then releasing the button to lock them in place (Fig. 32)

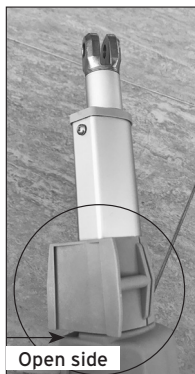
Fig. 33

8. Install the control unit on the actuator with the shortest cable. Perform the following operations before installing the control unit on the actuator:

- Connect the power cable to the control unit;
- Insert the plastic support onto the actuator, as shown in Fig. 33, with the open side of the tapered part facing downwards;
- Install the control unit on the actuator by sliding it into the guides and locking it in place using the plastic support installed in the previous point. Fig. 34;
- Install the actuator with the control unit on the two welded pins underneath the mesh structure at the headboard end. As follows:
  - Align the hole on the actuator with the one on the brackets of the bed section (Fig. 35);
  - Insert the quick-release pin into the hole (ensure the bracket is properly positioned inside the actuator) (Fig. 36);
  - Close the stop mechanism on the quick-release pin (Fig. 37).

Fig. 34

9. Assemble the other actuator underneath the structure of the mesh at the leg raiser end, repeating the operations described in the previous point.

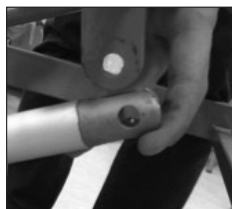
Fig. 35

Fig. 36

Fig. 37

10. Make the electrical connections on the control unit, as follows:

- Remove the plastic cover from the control unit and insert a flat screwdriver into the small window on the cover (Fig. 38);
- connect the cable of the actuator installed on the control unit on the "back" channel;
- connect the cable of the actuator on the foot raiser side on the "foot" channel;
- connect the cable from the actuator installed on the bed on the "Hi-Low" B channel;

Fig. 38

- A - Not connected
- B - Mesh lifting electric actuator
- C - Leg raiser electric actuator
- D - Head raiser electric actuator
- E - Remote control

- Position the cables in the cable housings present on the control unit and reposition the plastic cover

At this point, the bed is assembled and ready for use.

### 8.3 MP390 and MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA)

1. Insert the two rectangular tubes into the frame at the head raiser end and push them right in (Fig. 39)

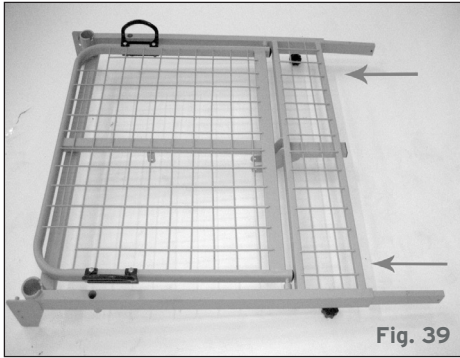
Fig. 39

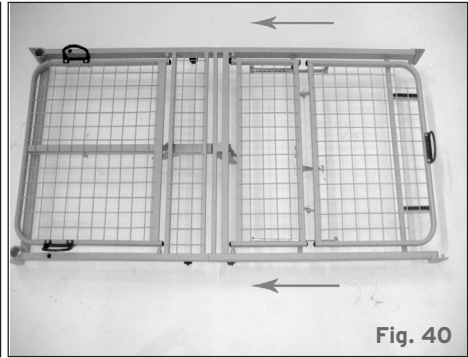
Fig. 40

Fig. 41

2. Insert the headboard frame and the footboard frame into the supports on the mesh frame, making sure that the locking pin is correctly inserted into the hole (Fig. 41). The headboard is the one with the shorter actuator cable.

3. Prepare the sides for assembly in the following manner:

- Take a pair of slats and position them parallel to the floor.

N.B. THE SIDES HAVE A SMOOTH SIDE, A LOWER PART AND A ROUNDED SIDE THAT IS THE UPPER PART. POSITION THE SMOOTH SIDE OF BOTH IN THE SAME DIRECTION.

Fig. 42

- Insert the plastic slides as a connection element of the sides on both sides, taking care to ensure that the pointed triangular part remains in the same direction as the rounded part of the side (Fig. 42)

- Carry out the same operation for the other two slats, in a manner to form the pair of sides

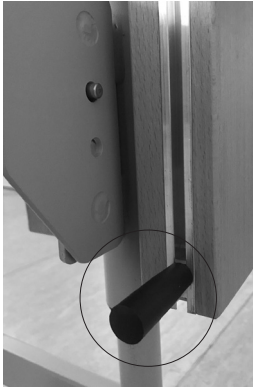
Fig. 43

4. Insert the two pairs of slats into the lower part of the guides on the wooden headboard and footboard. Perform this operation first on one side and then on the other. Tighten the stops at the end of the guides to prevent the sides from coming out of the guides (Fig. 43).

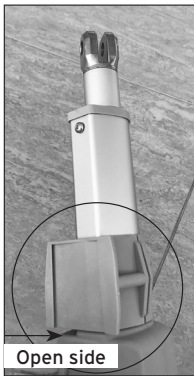
Fig. 44

5. Install the control unit on the actuator with the shortest cable. Perform the following operations before installing the control unit on the actuator:

- Connect the power cable to the control unit;
- Insert the plastic support onto the actuator, as shown in Fig. 44, with the open side of the tapered part facing downwards;
- Install the control unit on the actuator by sliding it into the guides and locking it in place using the plastic support installed in the previous point. Fig. 45;
- Install the actuator with the control unit on the two welded pins underneath the mesh structure at the headboard end. As follows:
  - Align the hole on the actuator with the one on the brackets of the bed section (Fig. 46);
  - Insert the quick-release pin into the hole (ensure the bracket is properly positioned inside the actuator) (Fig. 47);
  - Close the stop mechanism on the quick-release pin (Fig. 48).

Fig. 45

6. Assemble the other actuator underneath the structure of the mesh at the leg raiser end, repeating the operations described in the previous point.

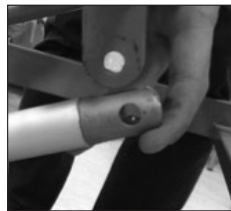
Fig. 46

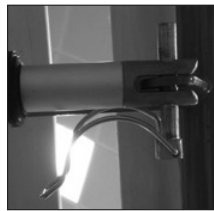
Fig. 47

Fig. 48

7. Make the electrical connections on the control unit, as follows:

- Remove the plastic cover from the control unit and insert a flat screwdriver into the small window on the cover (Fig. 49);

- connect the cable of the actuator installed on the control unit on the "back" channel;
- connect the cable of the actuator on the foot raiser side on the "foot" channel;
- connect the cable from the actuator installed on the bed on the "Hi-Low" B channel;

Fig. 49

- A - Electric actuator footboard end
- B - Electric actuator headboard end
- C - Leg raiser electric actuator
- D - Head raiser electric actuator
- E - Remote control

- Position the cables in the cable housings present on the control unit and reposition the plastic cover

At this point, the bed is assembled and ready for use.

#### 8.4 Assembly of the handles on the sections (for MP370/MP380/MP381/MP390/MP391)

1. Use the handle with screws/bolts provided for assembly of the handle on the bed sections (Fig. 50).
2. Assemble the handles using the holes present on the structure, at the points shown in the photo (Fig. 51)

Fig. 50

Illustrative images

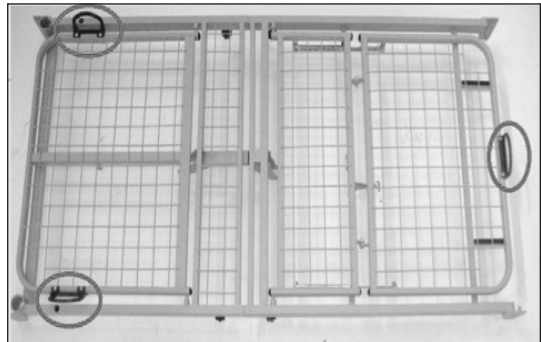
Fig. 51

### 9. BEFORE EVERY USE

- Always check for wear and tear of moving mechanical parts, in order to ensure that the wheelchair can be used in total safety for both people and property
- Check that the device is correctly assembled, taking particular care over the locking screws of the headboard, the footboard and the mechanical joints in general

**10. WARNINGS ON USE**

- Ensure that all the wheels are locked with their brake before making any adjustment to the bed
- Always disconnect the power cable from the mains electricity before performing any cleaning or maintenance operations or simply moving the bed
- Do not use additional mobile multiple plug sockets or extension cables
- Do not connect any electrical devices not intended for normal functioning to the bed

**WARNING!!!** The bed includes electrically controlled, moving parts. Operating the moving parts of the bed inappropriately could result in crushing of parts of the patient's limbs. Extreme care must be taken in using the electrical movements, making sure that the patient is positioned with all limbs stably resting on the top of the bed

**11. HOW TO USE**

**11.1 MP370 (three-jointed electric hospital bed mod. GERANIO)**


The MP370 (GERANIO) hospital bed offers the possibility of adjusting the inclination of the head raiser and the leg raiser, using the relative buttons on the remote control.

N.B. ALWAYS LOCK THE BED WHEELS WITH THE BRAKES ON EACH WHEEL BEFORE ANY ADJUSTMENT.

**WARNING**

Connect the power cable to a 220-230 VAC mains plug

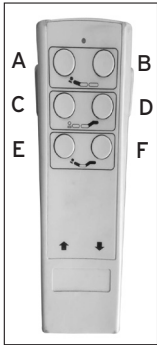- A Head raiser adjustment button (UP)
- B Head raiser adjustment button (DOWN)
- C Leg raiser adjustment button (UP)
- D Leg raiser adjustment button (DOWN)
- E Combined head raiser/leg raiser adjustment button (UP)
- F Combined head raiser/leg raiser adjustment button (DOWN)
- X Function blocking key (all bed functions are blocked if it is extracted)

N.B. THE DEVICE IN QUESTION MIGHT NOT BE INSTALLED ON THIS REMOTE CONTROL MODEL. THE MANUFACTURER RESERVES THE RIGHT TO DECIDE INDEPENDENTLY, BASED ON THE TYPE OF ELECTRONICS WITH WHICH THE BED IS EQUIPPED


**11.2 MP380-MP381 (three-jointed electric lifting hospital bed mod. GERBERA)**

The MP380 (GERBERA) hospital bed also offers the possibility of adjusting the inclination of the head raiser and the leg raiser, and also of raising the height of the mesh.

N.B. ALWAYS LOCK THE BED WHEELS WITH THE BRAKES ON EACH WHEEL BEFORE ANY ADJUSTMENT.

**WARNING**

Connect the power cable to a 220-230 VAC mains plug

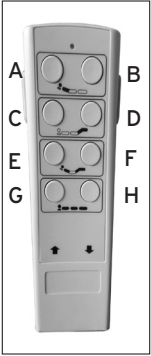- A Head raiser adjustment button (UP)
- B Head raiser adjustment button (DOWN)
- C Leg raiser adjustment button (UP)
- D Leg raiser adjustment button (DOWN)
- E Combined head raiser/leg raiser adjustment button (UP)
- F Combined head raiser/leg raiser adjustment button (DOWN)
- G Mesh lifting adjustment button (UP)
- H Mesh lifting adjustment button (DOWN)
- X Function blocking key (all bed functions are blocked if it is extracted)

N.B. THE DEVICE IN QUESTION MIGHT NOT BE INSTALLED ON THIS REMOTE CONTROL MODEL. THE MANUFACTURER RESERVES THE RIGHT TO DECIDE INDEPENDENTLY, BASED ON THE TYPE OF ELECTRONICS WITH WHICH THE BED IS EQUIPPED


### 11.3 MP390-MP391 (three-jointed electric Trendelenburg hospital bed mod. GARDENIA)

The MP390 (GARDENIA) hospital bed offers the possibility of adjusting the inclination of the head raiser and the leg raiser, raising the mesh and also adjusting the height of the headboard and footboard independently, thus allowing the bed to be placed in the Trendelenburg and reverse Trendelenburg position.

N.B. ALWAYS LOCK THE BED WHEELS WITH THE BRAKES ON EACH WHEEL BEFORE ANY ADJUSTMENT.

**WARNING**

Connect the power cable to a 220-230 VAC mains plug

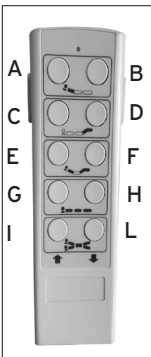- A Head raiser adjustment button (UP)
- B Head raiser adjustment button (DOWN)
- C Leg raiser adjustment button (UP)
- D Leg raiser adjustment button (DOWN)
- E Combined head raiser/leg raiser adjustment button (UP)
- F Combined head raiser/leg raiser adjustment button (DOWN)
- G Mesh lifting adjustment button (UP)
- H Mesh lifting adjustment button (DOWN)
- I "Reverse Trendelenburg" adjustment button (head UP, legs DOWN)
- L "Trendelenburg" adjustment button (head DOWN, legs UP)
- X Function blocking key (all bed functions are blocked if it is extracted)

N.B. THE DEVICE IN QUESTION MIGHT NOT BE INSTALLED ON THIS REMOTE CONTROL MODEL. THE MANUFACTURER RESERVES THE RIGHT TO DECIDE INDEPENDENTLY, BASED ON THE TYPE OF ELECTRONICS WITH WHICH THE BED IS EQUIPPED

## 12. MAINTENANCE

MOPEDIA devices manufactured by Moretti S.p.A. are CE marked and carefully checked before being commercialised. For patient and doctor safety, have the product checked at least once a year by the manufacturer or by an authorised laboratory to verify its suitability for use. Periodic inspection means examination of the following parts:

- The bearing structure of the bed and the lifting mechanism with relative attachments, controls, safety devices and patient support devices
- The sides of the bed (check attachment and movement)
- Check wheel attachment and brake functioning
- Check that the electric cables and plugs are undamaged
- Check all welded points

Use only original replacement parts and accessories for repairs.

## 13. CLEANING AND DISINFECTION

### 13.1 Cleaning

Use only a damp cloth and mild soap. Then dry thoroughly before use. Use water at a maximum temperature of 30 °C. Do not use water jet and/or steam cleaning machines.

### 13.2 Disinfection

Use a common disinfectant cleaner if you need to disinfect the wheelchair.

N.B. NEVER USE ACID, ALKALI OR SOLVENTS (ACETONE OR DILUENT) TO CLEAN IT.

## 14. CONDITIONS OF DISPOSAL

### 14.1 General conditions of disposal

Never dispose of the product as normal domestic waste. Dispose of the product at a sorted waste collection centre for recycling.

### 14.2 INSTRUCTIONS FOR CORRECT DISPOSAL IN ACCORDANCE WITH EUROPEAN DIRECTIVE 2002/96/EC:

At the end of its working life, the product must not be disposed of together with normal urban waste. It must be delivered to municipal separated waste collection facilities, or to appropriate dealers that provide this service. Separated waste disposal helps to reduce possible negative effects on the environment and health deriving from improper disposal and allows for recycling the materials comprising the product, which translates into significant energy and resource savings. The product bears the barred bin symbol to underline the obligation of disposing of electro-medical equipment.

## 15. DECLARATION OF ELECTROMAGNETIC COMPATIBILITY

The device is designed for use in the electromagnetic conditions specified in the following tables. The user of the device must make sure that these conditions obtain when he is using it.

Table 1

Manufacturer's guide and declaration - electromagnetic emissions		
The device is intended to be used in the following electromagnetic conditions. The client or user must ensure that it is used in these conditions.		
Emission tests	Conformity	Electromagnetic environment - guide
RF emissions CISPR 11	Group 1	The device uses RF energy only for internal functioning and for charging the battery. Its RF emissions are therefore very low and in all likelihood do not interfere with nearby electronic equipment.
Emissions RF CISPR 11	Class B	The device is intended to be used in all buildings, including the home and buildings connected directly to the low voltage public power grid for residential applications.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage flutter/flicker IEC 61000-3-3	Conforming	

Table 2

Manufacturer's guide and declaration - electromagnetic emissions			
The device is intended to be used in the following electromagnetic conditions. The client or user must ensure that it is used in these conditions.			
Immunity tests	Test level IEC 60601	Conformity level	Electromagnetic environment - guide
Electrostatic discharge (ESD) IEC 61000-4-2	+ 6 kV contact + 8 kV air	+ 6 kV contact + 8 kV air	Information for using the device and charging the battery. The floor must be wood, concrete or ceramic. If the floor is covered with a synthetic material, the relative humidity must be at least 30%.
Electrical fast transient/ burst immunity IEC 61000-4-4	+2KV for power lines +1 kV for I/O lines	2KV for power lines +1 kV for I/O lines	The quality of the mains voltage must be that of a typical business environment or hospital when charging the battery
Over voltages IEC 61000-4-5	1 kV differential mode +2 KV common mode	1 kV differential mode +2 KV common mode	The quality of the mains voltage must be that of a typical business environment or hospital when charging the battery
Voltage drop outs, brief interruptions and variations of voltage on the power supply lines. IEC 61000-4-11	<5% UT (>95% UT drop out) for 0.5 cycles 40% UT (60% UT drop out) for 5 cycles 70% UT (30% UT drop out) for 25 cycles <5% UT (>95% UT drop out) for 5 s.	<5% UT (>95% UT drop out) for 0.5 cycles 40% UT (60% UT drop out) for 5 cycles 70% UT (30% UT drop out) for 25 cycles <5% UT (>95% UT drop out) for 5 s.	The quality of the mains voltage must be that of a typical business environment or hospital when charging the battery.
Magnetic field Mains frequency (50/60 Hz) IEC 61000-4-8	3 A/m	3 A/m	Magnetic fields at the mains frequency should have the typical characteristics of a commercial or hospital environment. (Applies both to using the device and charging the battery).
Note: UT is the AC mains voltage prior to application of the test level.			

**Table 3**

Manufacturer's guide and declaration - electromagnetic immunity			
The device is intended to be used in the following electromagnetic conditions. The client or user must ensure that it is used in these conditions. Make sure it is used in these conditions.			
Immunity tests	Test level IEC 60601	Conformity level	Electromagnetic environment - guide
conducted RF IEC 61000-4-6	3 V eff. From 150 KHz to 80 MHz	3 V eff.	
radiated RF IEC 61000-4-3	3 V/m From 80 MHz to 2.5 GHz	3 V/m	<p>Portable and mobile RF communications devices should not be used near any part of the device, including the cables, unless the recommended separation distances are observed, calculated from the equation applicable to the transmitter frequency.</p> <p>Recommended distances d = 1.2√P</p> <p>d = 1.2√P From 80 MHz to 800 Mhz</p> <p>d = 2.3√P from 800 MHz to 2.5 GHz</p> <p>d = 7.7 m during charging operations</p> <p>where "P" indicates the maximum rated output of the transmitter in Watt (W) according to the transmitter manufacturer and "d" is the recommended separation distance in metres (m). The intensity of the field generated by fixed RF transmitters, as determined in an on-site electromagnetic inspection, might be below the level of conformity in each frequency band. Interference may be experienced in the proximity of equipment marked with this symbol:</p> 
<p>NOTE 1: At 80 MHz and 800 MHz, apply the separation distance for the higher frequency band.</p> <p>NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>			
<p><sup>a</sup> The intensity of the field for fixed transmitters, such as base stations for radio-telephones (cellulars and cordless) and terrestrial mobile radios, amateur radio equipment, AM and FM radio transmitters and TV transmitters may not be known precisely in theory. To establish an electromagnetic environment caused by fixed RF transmitters, one must run an on-site electromagnetic inspection. If the field intensity measured in the place where the DEVICE is used exceeds the applicable RF level of conformity indicated above, normal functioning of the DEVICE must be monitored. If abnormalities are observed, supplementary measures may be necessary, such as re-orientation or repositioning of the DEVICE.</p> <p><sup>b</sup> The field intensity over the frequency range from 150 kHz to 80 MHz should be less than 3 V/m.</p>			

Table 4

Recommended separation distances between portable and mobile radio communications equipment and the DEVICE						
<p>The device is intended to operate in an electromagnetic environment in which RF disturbances are under control. The client or user of the DEVICE can contribute to preventing EM interference by ensuring a minimum distance between the mobile and portable RF communications devices (transmitters) and the DEVICE battery charger during charging and the entire device during its operation as recommended below, in relation to the maximum output power of the radio communication equipment</p>						
Tests	Separation distance as a function of the frequency of the transmitter m					
Maximum rated output power of transmitter W	During battery charging	During use OF THE DEVICE	During battery charging	During use OF THE DEVICE	During battery charging	During use OF THE DEVICE
	from 150 kHz to 80 MHz $d = 1.2\sqrt{P}$	from 150 kHz to 80 MHz $d = 1.2\sqrt{P}$	from 80 kHz to 800 MHz $d = 1.2\sqrt{P}$	from 80 kHz to 800 MHz $d = 1.2\sqrt{P}$	from 800 kHz to 2.5 GHz $d = 1.2\sqrt{P}$	from 800 kHz to 2.5 GHz $d = 1.2\sqrt{P}$
0.01	0.12	0.12	0.12	0.12	0.23	0.23
0.1	0.38	0.38	0.38	0.38	0.73	0.73
1	1.2	1.2	1.2	1.2	2.3	2.3
10	3.8	3.8	3.8	3.8	7.3	7.3
100	12	12	12	12	23	23
<p>For transmitters whose maximum rated output power is not listed above, the recommended separation distance <math>d</math>, in metres (m), can be calculated using the equation applicable to the transmitter frequency, where <math>P</math> is the maximum rated output power of the transmitter, in Watt (W), as given by the transmitter's manufacturer.</p> <p>NOTE 1: At 80 MHz and 800 MHz, apply the separation distance for the higher frequency band.</p> <p>NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.</p>						

## 16. REPLACEMENT PARTS AND ACCESSORIES

Use only replacement parts and accessories listed in the Moretti general catalogue. Available accessories:

- MPA301 Patient lifter
- MPA303 Drip rod
- MPA305 Hospital bed table
- MPA307 Extension for sides

## 17. TROUBLESHOOTING

PROBLEMS	PROBABLE CAUSES	SOLUTION
<b>NO CONTROL IS FUNCTIONING</b>	<ol style="list-style-type: none"> <li>1. The electric plug is not connected to the mains</li> <li>2. The remote control is not connected to the control unit</li> <li>3. Probable fault on the remote control or control unit</li> </ol>	<ol style="list-style-type: none"> <li>1. Connect the electric plug</li> <li>2. Check and, if necessary, connect the remote control to the control unit</li> <li>3. Contact the assistance centre for a precise diagnosis of the problem</li> </ol>
<b>THE ELECTRIC HEAD RAISER IS NOT WORKING</b>	<ol style="list-style-type: none"> <li>1. The head raiser actuator is not connected to the control unit</li> <li>2. The actuator or the remote control or the control unit may be faulty</li> </ol>	<ol style="list-style-type: none"> <li>1. Check and, if necessary, connect the actuator to the control unit</li> <li>2. Contact the assistance centre for a precise diagnosis of the problem</li> </ol>
<b>THE ELECTRIC LEG RAISER IS NOT WORKING</b>	<ol style="list-style-type: none"> <li>1. The leg raiser actuator is not connected to the control unit</li> <li>2. The actuator or the remote control or the control unit may be faulty</li> </ol>	<ol style="list-style-type: none"> <li>1. Check and, if necessary, connect the actuator to the control unit</li> <li>2. Contact the assistance centre for a precise diagnosis of the problem</li> </ol>
<b>ELECTRIC LIFTING OF THE BED IS NOT WORKING (ONLY MP380 AND MP381)</b>	<ol style="list-style-type: none"> <li>1. The lifting actuator is not connected to the control unit</li> <li>2. The actuator or the remote control or the control unit may be faulty</li> </ol>	<ol style="list-style-type: none"> <li>1. Check and, if necessary, connect the actuator to the control unit</li> <li>2. Contact the assistance centre for a more precise diagnosis of the problem</li> </ol>
<b>TRENDELENBURG/ REVERSE TRENDELENBURG LIFTING OF THE BED IS NOT WORKING (ONLY MP390 AND MP391)</b>	<ol style="list-style-type: none"> <li>1. One or both of the actuators (headboard/footboard) are not connected to the control unit</li> <li>2. The actuator(s) or the remote control or the control unit may be faulty</li> </ol>	<ol style="list-style-type: none"> <li>1. Check and, if necessary, connect the actuators to the control unit</li> <li>2. Contact the assistance centre for a more precise diagnosis of the problem</li> </ol>

## 18. TECHNICAL SPECIFICATIONS

### 18.1 Technical specifications for all models

<b>UNITS OF FORCE NEEDED TO CONTROL THE COMMANDS ON THE REMOTE CONTROL</b>	5N
<b>POWER SUPPLY</b>	220 VAC 50-60 Hz
<b>VOLTAGE OUTPUT FROM THE CONTROL UNIT</b>	24 VDC
<b>PROTECTION CLASS</b>	IP 54
<b>NOISE LEVEL</b>	< 45 dB
<b>INSULATION CLASS</b>	Class 2
<b>MATERIALS USED FOR MP370-380-381-390-391</b>	Powder-coated steel tube structure, welded steel mesh (MP370-380-390) or wooden base (MP381-391). Wooden headboard, footboard and sides.

### 18.2 Dimensions

#### MP370 Three-jointed electric hospital - bed (mod. GERANIO)

Bed dimensions: 208x102x94h cm (at the headboard and footboard)

Mesh dimensions: 187x80x42.5h cm (mesh base height from the floor)

Maximum height at the headboard/footboard: 94 cm

Maximum angle of the head raiser: 70°

Maximum angle of the leg raiser: 30°

<b>MAXIMUM WEIGHT</b>	135 Kg
<b>TOTAL WEIGHT</b>	77.5 Kg
<b>HEADBOARD WEIGHT</b>	19 Kg
<b>FOOTBOARD WEIGHT</b>	19 Kg
<b>FRAME WEIGHT HEAD RAISER END</b>	18 Kg
<b>FRAME WEIGHT LEG RAISER END</b>	16 Kg
<b>ELECTRICAL KIT WEIGHT</b>	5.5 Kg

**MP380 - MP381 Three-jointed electric lifting hospital bed (mod. GERBERA)**

Bed dimensions: 208x102x116.5h cm (at the headboard and footboard)

Mesh dimensions: 187x80x33.5+65.5h cm (mesh base height from the floor)

Maximum height at the headboard/footboard: 116.5 cm

Minimum height at the headboard/footboard: 22.5 cm

Minimum height of frame from the ground: 22.5 cm

Maximum angle of the head raiser: 70°

Maximum angle of the leg raiser: 30°

<b>MAXIMUM WEIGHT</b>	135 Kg
<b>TOTAL WEIGHT</b>	99 Kg
<b>HEADBOARD WEIGHT</b>	19 Kg
<b>FOOTBOARD WEIGHT</b>	19 Kg
<b>FRAME WEIGHT HEAD RAISER END</b>	18 Kg
<b>FRAME WEIGHT LEG RAISER END</b>	16 Kg
<b>FRAME WEIGHT ON LIFTING</b>	20 Kg
<b>ELECTRICAL KIT WEIGHT</b>	7 Kg

**MP390-MP391 Three-jointed electric Trendelenburg hospital bed (mod. GARDENIA)**

Bed dimensions: 216x102x127h cm

Mesh dimensions: 187x80x40 cm+81h

Maximum height at the headboard/footboard: 127 cm

Minimum height at the headboard/footboard: 25 cm

Maximum angle of the head raiser: 70°

Maximum angle of the leg raiser: 30°

Maximum "Trendelenburg" angle (head down, legs up): 12°

Maximum "reverse Trendelenburg" angle (head up, legs down): 12°

<b>MAXIMUM WEIGHT</b>	135 Kg
<b>TOTAL WEIGHT</b>	84 Kg
<b>HEADBOARD WEIGHT</b>	19 Kg
<b>FOOTBOARD WEIGHT</b>	19 Kg
<b>FRAME WEIGHT HEADBOARD</b>	18 Kg
<b>FRAME WEIGHT FOOTBOARD</b>	16 Kg
<b>ELECTRICAL KIT WEIGHT</b>	12 Kg

## 19. WARRANTY

Moretti products are guaranteed for 2 (two) years from the date of sale against material and manufacturing defects, subject the following limitations. The warranty is voided by improper use, abuse, modifications to the product and failure to follow the instructions. The intended use of the product is given in the user manual

Moretti is not liable for damage, injury or any other consequences resulting from installation or use which are not scrupulously conforming with the instructions given in the installation, assembly and user manual.

Moretti does not guarantee its products against damage or defects in the following circumstances: natural disasters, unauthorised repair or maintenance, improper electric power supply (as applicable), use of parts or components not supplied by Moretti, failure to follow the guidelines and instructions for use, tampering, shipping damage (other than the original shipping by Moretti), or failure to perform maintenance as indicated in the manual.

Components subject to wear and tear are not covered by this warranty if the damage is caused by normal use of the product.

## 20. REPAIRS

Repairs under warranty

If a Moretti product has material or manufacturing defects during the warranty period, Moretti will agree with the client whether the defect is covered by the warranty. Moretti, at its sole discretion, may replace or repair the article at a specified Moretti reseller or its own premises. The costs of labour incurred in repairing the product will be borne by Moretti if it determines that the repair is covered by the warranty. Repair and replacement do not renew the warranty period.

Repairing a product not covered by the warranty

A product not covered by warranty may be returned for repair only if authorised in advance by Moretti customer service. The costs of labour and shipping incurred by repairs not covered by the warranty are borne by the client or reseller in their entirety. Repairs on products not covered by the warranty are themselves guaranteed for 6 (six) months from the day of reception of the repaired product.

Non-defective products

The client will be notified if Moretti concludes that the product is not defective after having received and examined it. The product will be returned to the client at his expense.

## REPLACEMENT PARTS

Moretti original replacement parts are guaranteed for 6 (six) months from the day of delivery.

## 21. NON-LIABILITY CLAUSE

Unless otherwise expressly specified in this warranty and within the limits of the law, Moretti makes no declaration, guarantee or condition, express or implicit, including any future declaration, guarantee or condition of sale, suitability for a given purpose, non violation and non interference. Moretti does not guarantee that the use of its product will be uninterrupted and problem-free. The duration of any implicit guarantee under the law is limited to the warranty period, within the limits of the law. Certain states and countries do not permit limitations on the duration of an implicit guarantee or the exclusion of limitation of accidental or indirect damages in relation to consumer products. In said states and countries, certain exclusions and limitations of this warranty may not apply to the user. This warranty is subject to modification without notice.

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